


## Configuration of SEEONE VIEW

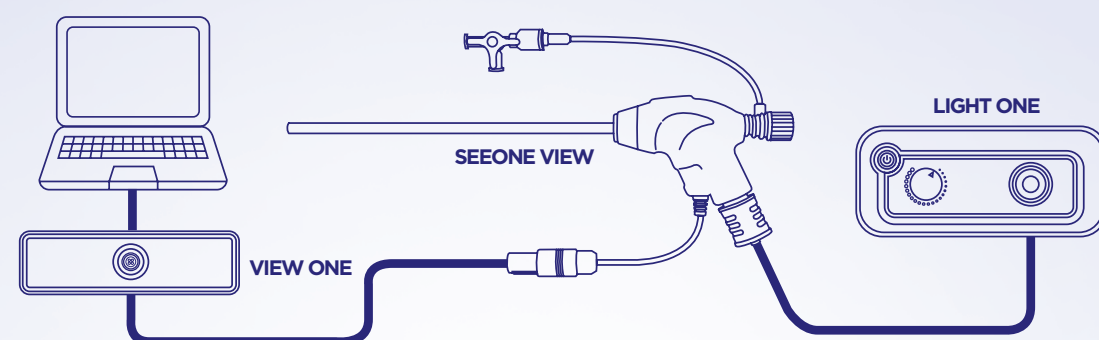

| View One Specification |         |                         | Light One Specification |             |                    |
|------------------------|---------|-------------------------|-------------------------|-------------|--------------------|
| Type                   |         | Image processing device | Type                    |             | Endolight          |
| Size                   |         | 150X140X45              | Size                    |             | 210X255X110        |
| Resolution             |         | 160k (400x400 pixel)    | Power Supply            | Voltage     | 100-240VAC         |
| Power Supply           | Voltage | USB 5Vdc                |                         | Frequency   | 50/60Hz            |
|                        | Power   |                         | LED                     | Power       | 60W (15A)          |
| Cable                  | Video   | 2m                      |                         | CCT range   | 6500k (Cold white) |
|                        | Type    | USB C 2.0 compatible    |                         | Brightness  | 0~100              |
|                        |         |                         | Cable                   | Fiber optic | 2m                 |
|                        |         |                         |                         | Power       | 1.5m               |

# SEEONE VIEW

## Epidural Endoscopic Catheter System


- Microendoscopic camera with high resolution & 160K
- Accurate treatment with visualization
- Minimally invasive

SEEONE VIEW

Epidural Endoscopic Catheter System

SEEONE VIEW is a micro-endoscope system with 160K high-definition camera, working channel and optical fiber for LED light.

It provides real-time images using the high-definition camera by delivering light supplied from light source through optical fiber. Therefore, the largest benefit of SEEONE VIEW is that users can directly check the region of lesion and make a precise diagnosis and treatment. Moreover, drugs can be injected through the working channel.

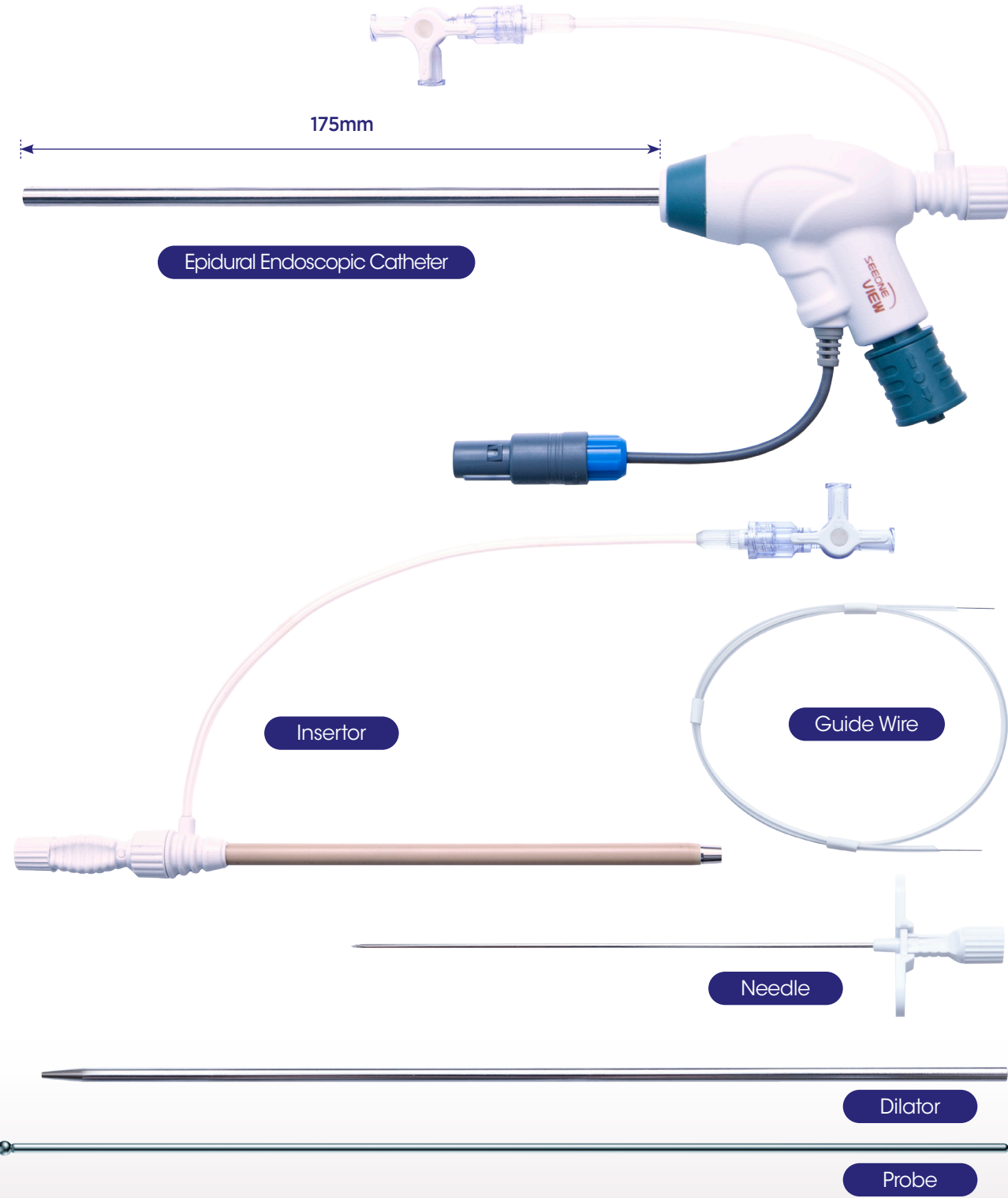

Features



High-Resolution

High Resolution Camera

The high- resolution enables precise diagnosis and surgical operation, with more rapid access to lesions and reducing surgery time

Safe

Minimally invasive

Minimum incision to reduce pain and recovery time.

Accurate treatment with visualization and illumination

Provides superior images enabling more safe operation.

Easy

Easy-to-use program


It uses a dedicated program

Easy-to-use ergonomic design


Ergonomic handle design easy to use

View device for illumination intensity control

Convenient for users to directly adjust illumination intensity on the view.

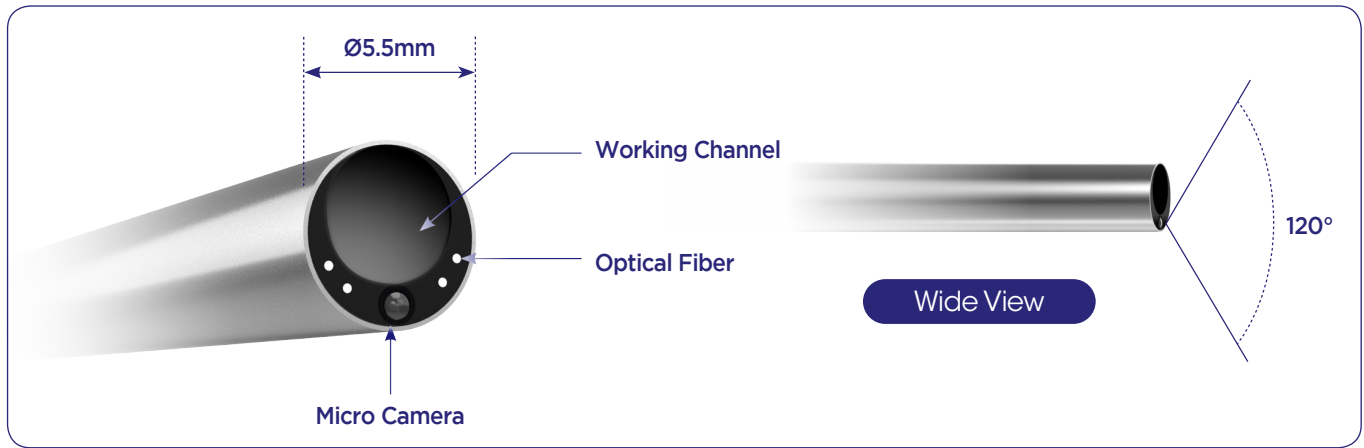

Optimal image quality

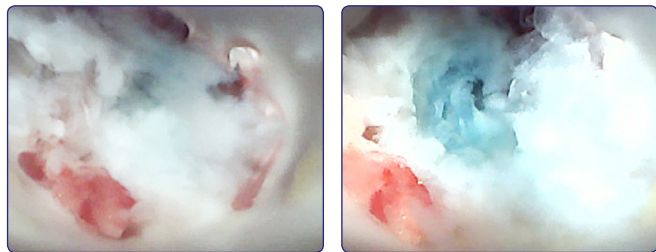

Model Spec

| Epidural Endoscopic Catheter |                     |
|------------------------------|---------------------|
| Type                         | Endoscope           |
| Endoscope Diameter           | Ø5.5mm              |
| Working Channel Diameter     | Ø3.0mm              |
| Angle of View                | 120°                |
| Angle of Curve               | 0°                  |
| Effective Pixels             | 160,000 (400 x 400) |